

ATW0036

Aptamer to Bovine IFN γ

Selection Information

Target for Selection: Bovine interferon gamma, GenWay Biotech # GWB-0F0115

Number of DNA Nucleotides: 32

Aptamers were selected from a randomized Base Pair 32-mer DNA library against the target protein. Proprietary methods were used to select this specific aptamer sequence.

Affinity Determination

Affinity Determination Method: Microscale Thermophoresis (MST)

Buffer Used for Affinity Determination: 1 x PBS, 1 mM MgCl₂, pH 7.4

Average kD: 33.7 \pm 6.0 nM

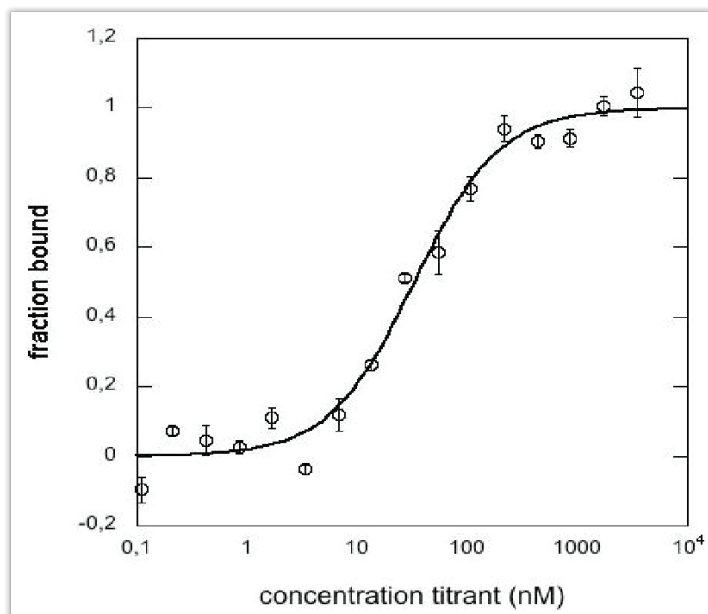


Figure 1. Aptamer-IFN γ Binding

The fraction of bound aptamer is plotted versus the titrated target concentration, $R^2 = 0.981$

Aptamer Folding

For optimal binding, aptamers must be folded into their tertiary structure prior to use. Once the aptamer is at its working concentration, heat to 90-95°C for 5 minutes, then cool to room temperature (~15 minutes)

For Research Use Only