

ATW0035

Aptamer to Human IL-6

Selection Information

Target for Selection: Recombinant Human IL-6, R&D Systems Cat# 206-IL-010

Number of DNA Nucleotides: 32

Aptamers were selected from a randomized Base Pair 32-mer DNA library against the target protein. Proprietary methods were used to select this specific aptamer sequence.

Affinity Determination

Affinity Determination Method: Microscale Thermophoresis analysis (MST)

Buffer Used for Affinity Determination: 1 x PBS, 1mM MgCl₂, 0.05% Tween 20, pH 7.4

Average kD: 19 ± 4.5 nM

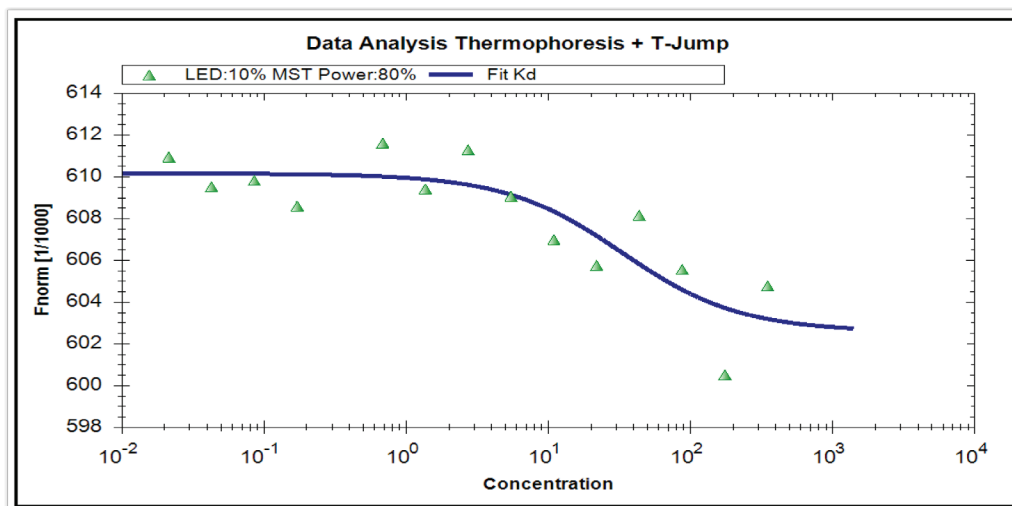


Figure 1. Aptamer-IL-6 Binding

The thermophoresis is plotted versus the titrated IL-6 concentration

Aptamer Folding

For optimal binding, aptamers must be folded into their tertiary structure prior to use. Once the aptamer is at its working concentration, heat to 90-95°C for 5 minutes, then cool to room temperature (~15 minutes)

For Research Use Only