

## ATW0052

Aptamer to Human CEACAM-1 / CD66a

### Selection Information

**Target for Selection:** Human CEACAM-1 / CD66a Protein, R&D Systems Cat# 2244-CM-050

**Number of DNA Nucleotides:** 32


Aptamers were selected from a randomized Base Pair 32-mer DNA library against the target protein. Proprietary methods were used to select this specific aptamer sequence.

### Affinity Determination

**Affinity Determination Method:** Microscale Thermophoresis analysis (MST)

**Buffer Used for Affinity Determination:** 1 x PBS, 1mM MgCl<sub>2</sub>, 0.05% Tween 20, pH 7.4

**Average kD:** 98.4 ± 15.3 nM

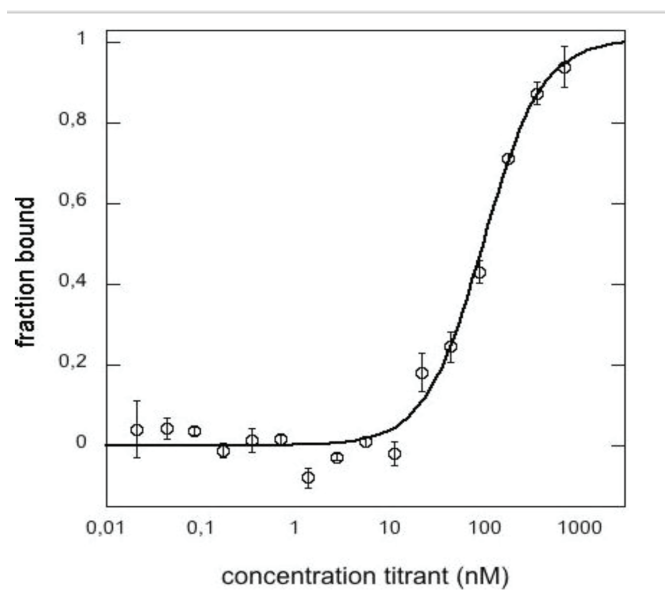

**Figure 1. Aptamer-CD66a Binding**

The fraction of bound aptamer is plotted versus the titrated target concentration,  $R^2 = 0.986$

### Aptamer Folding

For optimal binding, aptamers must be folded into their tertiary structure prior to use. Once the aptamer is at its working concentration, heat to 90-95°C for 5 minutes, then cool to room temperature (~15 minutes)

*For Research Use Only*