

ATW0053

Aptamer to Cytokeratin 19

Selection Information

Target for Selection: Recombinant human Cytokeratin 19 (KRT19), ProSopec Bio Cat# pro-350

Number of DNA Nucleotides: 70 (with primers)

Aptamers were selected from a randomized Base Pair 32-mer DNA library against the target protein. Proprietary methods were used to select this specific aptamer sequence.

Affinity Determination

Affinity Determination Method: Microscale Thermophoresis analysis (MST)

Buffer Used for Affinity Determination: 1 x PBS, 1mM MgCl₂, 0.05% Tween 20, pH 7.4

Average kD: 579.1 ± 138.6 nM

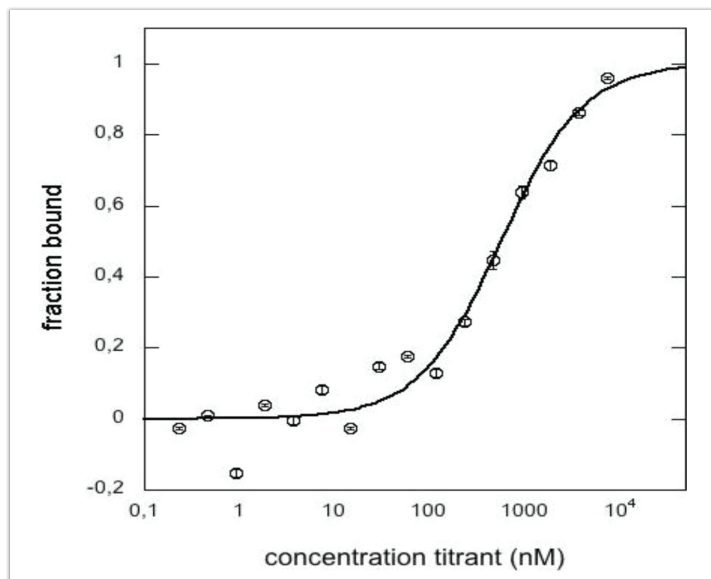


Figure 1. Aptamer-Cytokeratin 19 Binding

The fraction of bound aptamer is plotted versus the titrated target concentration, $R^2 = 0.97$

Aptamer Folding

For optimal binding, aptamers must be folded into their tertiary structure prior to use. Once the aptamer is at its working concentration, heat to 90-95°C for 5 minutes, then cool to room temperature (~15 minutes)

For Research Use Only